


Mix and Dispenche Protocols

Liquidhandling Robot Deck

Liquidhandling Robot Deck:

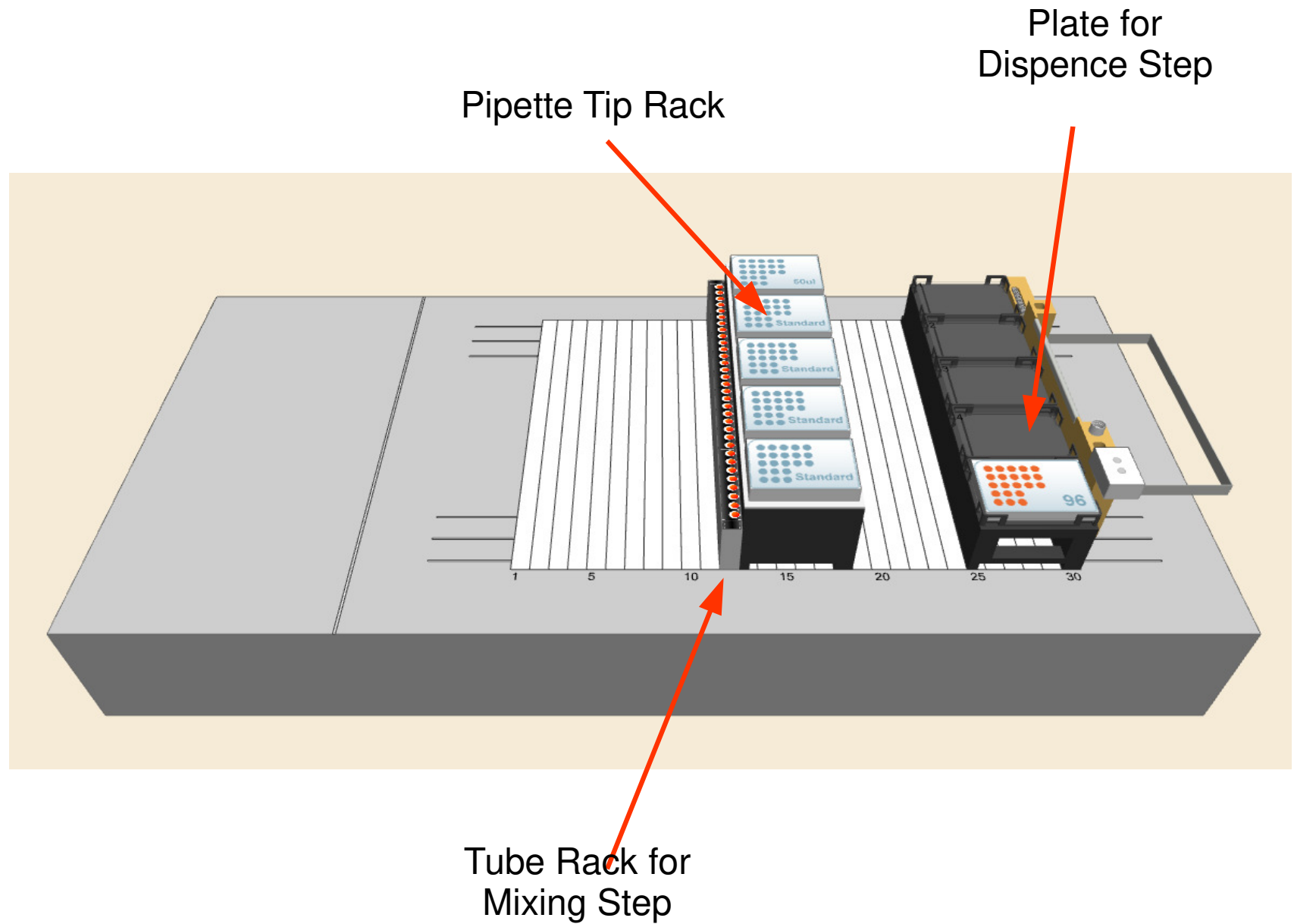

Recipe File (xlsx)

Recipe File for Mixing:

Workbook Name: Mix

Source Tube
Carrier
1. Component

Source Tube
Position
1. Component

Mixing Ratio

Final Volume

A	B	C	D	E	F	G	H	I	J
Command	ChannelPattern	SourceLabwareID_1	SourcePosID_1	SourceRatio_1	SourceLabwareID_2	SourcePosID_2	TargetLabwareID	TargetPosID	Volume
Chg	1111	SampleCarrier_Eppendorf	28	0.5	SampleCarrier_Eppendorf	27	SampleCarrier_Eppendorf	22	600
Chg	1111	SampleCarrier_Eppendorf	29	0.5	SampleCarrier_Eppendorf	26	SampleCarrier_Eppendorf	21	600
Chg	1111	SampleCarrier_Eppendorf	30	0.5	SampleCarrier_Eppendorf	29	SampleCarrier_Eppendorf	20	600
Chg	1111	SampleCarrier_Eppendorf	31	0.5	SampleCarrier_Eppendorf	30	SampleCarrier_Eppendorf	19	600
Chg	1111	SampleCarrier_Eppendorf	28	0.5	SampleCarrier_Eppendorf	25	SampleCarrier_Eppendorf	18	600
Chg	1111	SampleCarrier_Eppendorf	29	0.5	SampleCarrier_Eppendorf	28	SampleCarrier_Eppendorf	17	600
Chg	1111	SampleCarrier_Eppendorf	30	0.5	SampleCarrier_Eppendorf	27	SampleCarrier_Eppendorf	16	600
Chg	1111	SampleCarrier_Eppendorf	31	0.5	SampleCarrier_Eppendorf	28	SampleCarrier_Eppendorf	15	600
Chg	1111	SampleCarrier_Eppendorf	27	0.5	SampleCarrier_Eppendorf	26	SampleCarrier_Eppendorf	14	600
Chg	1111	SampleCarrier_Eppendorf	29	0.5	SampleCarrier_Eppendorf	24	SampleCarrier_Eppendorf	13	600
Chg	1111	SampleCarrier_Eppendorf	30	0.5	SampleCarrier_Eppendorf	25	SampleCarrier_Eppendorf	12	600
Chg	1111	SampleCarrier_Eppendorf	31	0.5	SampleCarrier_Eppendorf	26	SampleCarrier_Eppendorf	11	600
Chg	1111	SampleCarrier_Eppendorf	27	0.5	SampleCarrier_Eppendorf	24	SampleCarrier_Eppendorf	10	600
Chg	1111	SampleCarrier_Eppendorf	26	0.5	SampleCarrier_Eppendorf	25	SampleCarrier_Eppendorf	9	600
Chg	1111	SampleCarrier_Eppendorf	25	0.5	SampleCarrier_Eppendorf	24	SampleCarrier_Eppendorf	8	600
Chg	1111	SampleCarrier_Eppendorf	31	0.5	SampleCarrier_Eppendorf	24	SampleCarrier_Eppendorf	7	600

Disable Tip
Exchange
FOR TESTING
ONLY

Source Tube
Carrier
2. Component

Source Tube
Position
2. Component

Target Tube
Carrier

Target Tube
Position

Recipe File for Dispence Step:

Workbook Name: Dispence

Volume

Source Tube
Carrier

Source Tube
Position


Target Tube
Carrier

Target Tube
Position

D	E	F	G	H
SourceLabwareID	SourcePosID	TargetLabwareID	TargetPosID	Volume
SampleCarrier_Eppendorf	31	Target_1	A1	100
SampleCarrier_Eppendorf	30	Target_1	A2	100
SampleCarrier_Eppendorf	29	Target_1	A3	100
SampleCarrier_Eppendorf	28	Target_1	A4	100
SampleCarrier_Eppendorf	27	Target_1	A5	100
SampleCarrier_Eppendorf	26	Target_1	A6	100
SampleCarrier_Eppendorf	25	Target_1	A7	100
SampleCarrier_Eppendorf	24	Target_1	A8	100
SampleCarrier_Eppendorf	19	Target_1	A9	100
SampleCarrier_Eppendorf	17	Target_1	A10	100
SampleCarrier_Eppendorf	14	Target_1	A11	100
SampleCarrier_Eppendorf	8	Target_1	A12	100
SampleCarrier_Eppendorf	26	Target_1	B1	100
SampleCarrier_Eppendorf	15	Target_1	B2	100
SampleCarrier_Eppendorf	21	Target_1	B3	100
SampleCarrier_Eppendorf	31	Target_1	B4	100
SampleCarrier_Eppendorf	25	Target_1	B5	100
SampleCarrier_Eppendorf	21	Target_1	B6	100
SampleCarrier_Eppendorf	15	Target_1	B7	100
SampleCarrier_Eppendorf	28	Target_1	B8	100
SampleCarrier_Eppendorf	17	Target_1	B9	100
SampleCarrier_Eppendorf	18	Target_1	B10	100

GUI Overview

Protocol Selecton (Main Dialog)

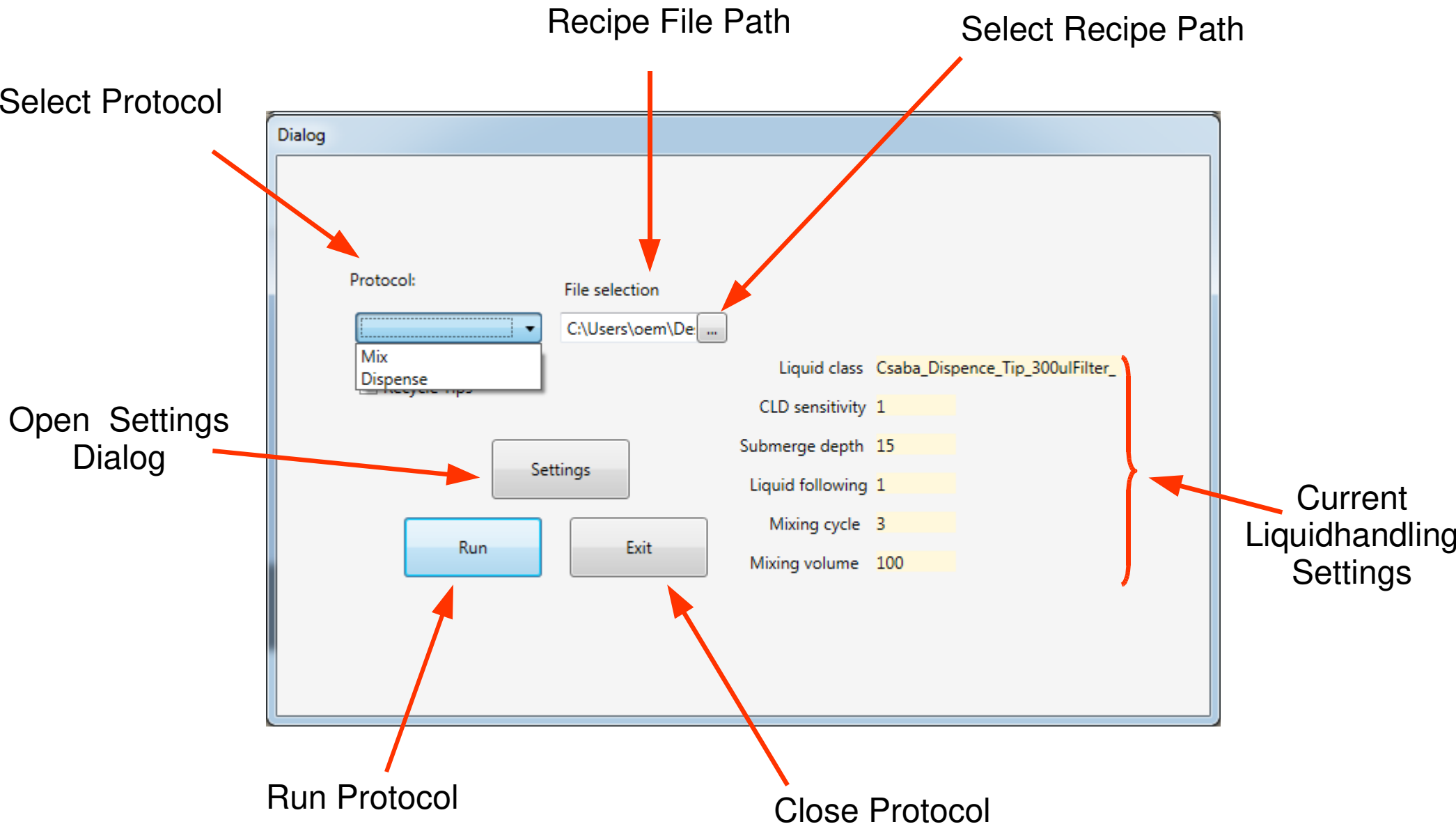

Liquid Handling Settings (Settings Dialog)

Dialog

Settings

Default

Liquid class
Csaba_Dispence_Tip_300

CLD Sensitivity
1

Submerge depth
15

Liquid following
1

Mixing cycle
3

Mixing volume
100

Save Exit